


Eliminate™ components

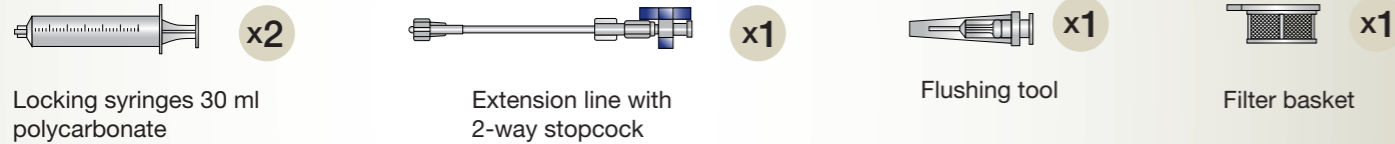

Ordering information

Guide catheter compatibility	Eliminate™
6Fr (I.D. ≥ 0.70 inch / 1.78 mm)	EG1602
7Fr (I.D. ≥ 0.80 inch / 2.03 mm)	EG1652

Product specifications

Tip Length	6 mm	Positioning marks	10 cm = single white mark located at 90 cm from distal tip
Coating	Hydrophilic	Wire compatibility	0.014" (0.36 mm)
Rapid exchange segment	23 cm	Catheter length	140 cm
Radiopaque marker band	1 mm located at 4 mm from the tip		



Terumo Corporation
44-1, 2-chome, Hatagaya
Shibuya-ku, Tokyo 151-0072,
Japan
Tel.: +81 3 3374 8111
Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 38 12 11
Fax: +32 16 40 02 49

**Terumo Europe
Interventional Systems
Liaison Office**
Immeuble Krebs
8, rue Jean-Pierre Timbaud
CS 10720
78067 Saint Quentin
en Yvelines Cedex,
France
Tel.: +33 139 30 95 75
Fax: +33 139 30 95 70

European Offices

**Terumo Europe N.V.
Benelux Sales Division**
Interleuvenlaan 40
3001 Leuven,
Belgium
Tel.: +32 16 39 25 80
Fax: +32 16 39 25 99
The Netherlands:
Tel.: 0800 0220396
Fax: 0800 0220414

Laboratoires Terumo France S.A.
Parc Ariane - Bât. Uranus
1, rue Hélène Boucher
78284 Guyancourt Cedex,
France
Tel.: +33 130 96 13 00
Fax: +33 130 43 60 85

Terumo Italia Srl
Via Simone Martini, 143 / 145
00142 Roma RM,
Italy
Tel.: +39 0651 96 14 20
Fax: +39 0650 30 407

Terumo UK Ltd
Tamesis, The Causeway
Egham, Surrey
TW20 9AW, United Kingdom
Tel.: +44 1784 476 200
Fax: +44 1784 476 201

Terumo Deutschland GmbH
Ludwig-Erhard-Strasse 6
65760 Eschborn, Germany
Tel.: +49 6196 80 230
Fax: +49 6196 80 23 200

Terumo Deutschland GmbH
Zweigniederlassung Spreitenbach
Bodenackerstrasse 3
8957 Spreitenbach, Switzerland
Tel.: +41 56 419 10 10
Fax: +41 56 419 10 11

Terumo Europe España SL
Avda. Juan Carlos I, Nº 13 - 7ª Planta
Edificio Torre La Garena
28806 Alcalá de Henares (Madrid),
Spain
Tel.: +34 9021 01 298
Fax: +34 9021 01 358

Terumo Sweden AB
Nya Varvet, Byggnad 90
Sven Källfelts Gata 18
426 71 Västra Frölunda, Sweden
Tel.: +46 3174 85 880
Fax: +46 3174 85 890

**Terumo Denmark
Filial of Terumo Sweden AB**
Tel.: +45 7020 93 80
Fax: +45 7020 94 80

**Terumo Europe N.V.
Moscow Representative Office**
JSC Northern Tower
13Fr., 10 Testovskaya Street
Moscow 123317, Russia
Tel.: + 7 495 988 4740
Fax: + 7 495 988 4739

**Addresses subject to change,
please consult our Website.**

© Registered Trademark
Published by Terumo Europe N.V.

Aspiration catheter

The deliverable solution for
thrombus removal

IS200GBE1109MDI03(071)JE

The deliverable solution for thrombus removal

Eliminate™ is an innovative aspiration catheter designed to offer an **optimal balance** between **crossing** performance, kink resistance and **thrombus aspiration capability**

Deliverability


Eliminate's innovative technology provides improved deliverability and kink resistance

- **Pre-loaded stylet** provides improved pushability and ensures optimal crossability in difficult to cross scenarios
- **Fully braided shaft construction** for improved kink resistance of the catheter
- **23 cm long rapid exchange segment** provides improved trackability
- **Hydrophilic coating** facilitates easy navigation through tortuous anatomies

Stiffening stylet with floppy distal segment



Fully braided shaft

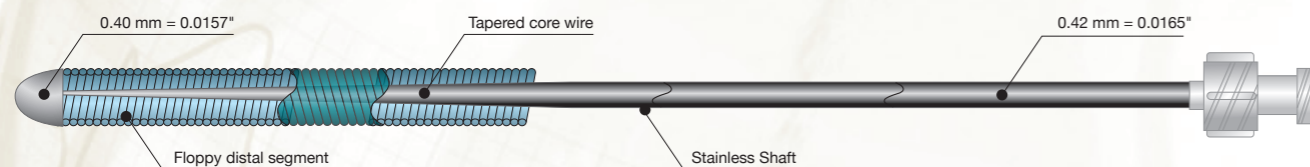

143 cm stiffening stylet


6 Fr guide catheter compatibility

Distal part



Extraction area
0.78 mm²


Proximal part

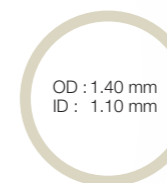

Extraction area
0.95 mm²


7 Fr guide catheter compatibility

Distal part

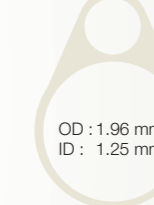

Extraction area
1.23 mm²

Proximal part

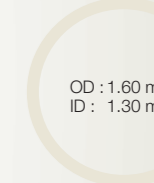


Extraction area
1.32 mm²

Aspiration Capability

Advanced technology for enhanced aspiration

- **Large extraction lumen** and the braided shaft construction ensures constant aspiration performance throughout the procedure
- **Rounded and short tip design** preserves crossability, while providing safe and effective aspiration
- A choice of **2 sizes** (6 Fr or 7 Fr guide catheter compatibility) to meet different clinical scenarios

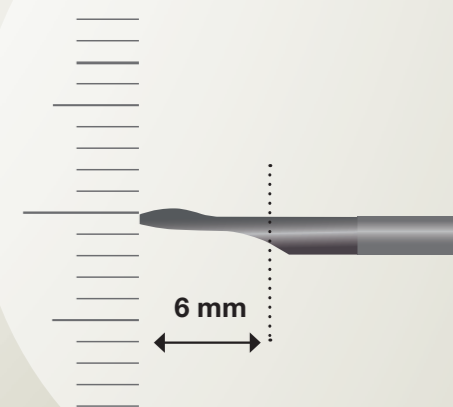

Ease of Use

Single-handed operator

- **Short and thin distal radiopaque marker** provides excellent distal visualization and increases distal flexibility
- **Complete kit contains:**
 - 30 ml locking syringes
 - Extension line (25 cm) with 2 way stopcock
 - Filter basket