

Technical specifications

Balloon material	Polyn P®
Balloon folding	3 and 6 folds*
Coating	ProGlide®
Balloon length	from 10mm to 200mm*
Balloon diameter	from 1.5 to 12mm*
Catheter length	80,120, 130 and 150cm*
Guidewire compatibility	.014" - .018" - .035"
Delivery Systems	.014" - Rx & OTW
	.018" OTW
	.035" OTW
Shelf life	3 years

Comparative analysis upon request

Contact us on info@raymed.ru for more information

Sheath & Guiding Catheter Compatibility

∅	2	3	4	5	6	7	8	9	10	12
Sheath	.035"			5F				6F		7F
	.018"		4F		5F					
	.014"		4F		5F					
Guiding	.035"			7F				8F		9F
	.018"		5F	6F	7F					
	.014"		4F	6F	7F					

Manufactured by Natec Medical Ltd

Distributed by

 Raymed Trading Group


129626, Moscow,
Grafskiy lane, 12a, b.2.

e-mail: info@raymed.ru
www.raymed.ru
+7 (495) 956 83 16



 **Raymed Trading Group**

Agile performance

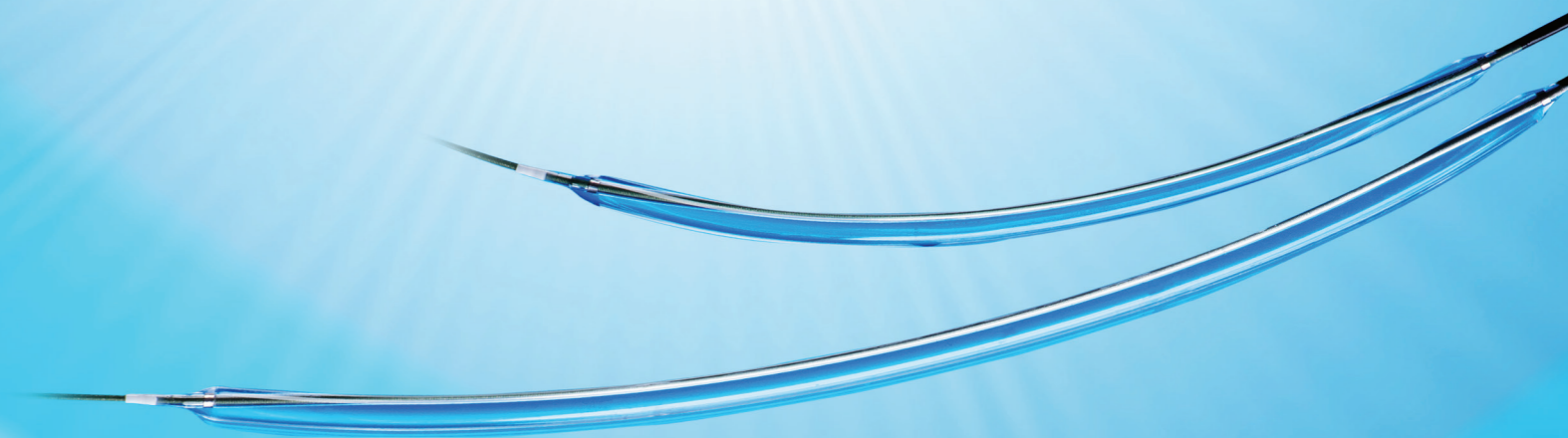

Enhanced crossability

Optimum track and push

Balanced performance

NIMBLE 

.035" OTW | .018" OTW | .014" OTW/RX

Nimble .035" OTW

Optimizing the performance by improving the deliverability

Because tip design makes the difference

- > Better flexibility along the guidewire
- > Perfect tip adherence on the guidewire (no fishmouth effect)
- > No kinking and no tip wrapping

Ordering information

Balloon Length (mm)	Balloon Diameter (mm)								
	3.00	4.00	5.00	6.00	7.00	8.00	9.00	10.00	12.00
20									
40									
60									
80									
100									
120									
150									
200									

Available sizes

Nimble .018" OTW

Reduce the friction... again again ... again ... again ...

- > Proprietary ProGlide® coating improving the trackability by 25% and reducing friction forces by 40% ... 4 times more durable than competitors**.

Unique tripod shaft construction

- > Reinforced shaft construction enabling more axial strength and higher pushability of the system
- > Comparable inflation and deflation time

Ordering information

Balloon Length (mm)	Balloon Diameter (mm)							
	2.00	2.50	3.00	3.50	4.00	5.00	6.00	7.00
40								
60								
80								
100								
120								
150								
200								

Available sizes

Nimble .014" OTW & Rx

- > Exceptional trackability on the longest sizes when balloon length can be longer than the shaft length outside the guiding catheter.

Infrapopliteal lesions and much more :

- > OTW & Rx shaft construction
- > 1.5 to 5.0mm balloons diameter
- > From 5.5 to 8.0 upon request for specific applications

Ordering information

Balloon Length (mm)	Balloon Diameter (mm)					
	1.50	2.00	2.50	3.00	3.50	4.00
20						
40						
60						
80						
100						
120						
150						
200						

Available sizes